



**Test and Evaluation Report  
of the  
Corometrics Medical Systems Neonatal  
Monitor Model 506**

**By**

**James E. Bruckart (Project Officer)  
Martin Quattlebaum (Project Officer)  
Brendan Squire (Project Officer)**

**Biodynamics Research Division**

**Bill Olding (UES)**

**August 1992**

**Approved for public release; distribution unlimited**

**United States Army Aeromedical Research Laboratory  
Fort Rucker, Alabama 36362-5292**

## Notice

### Qualified requesters

Qualified requesters may obtain copies from the Defense Technical Information Center (DTIC), Cameron Station, Alexandria, Virginia 22314. Orders will be expedited if placed through the librarian or other person designated to request documents from DTIC.

### Change of address

Organizations receiving reports from the U.S. Army Aeromedical Research Laboratory on automatic mailing lists should confirm correct address when corresponding about laboratory reports.

### Disposition

Destroy this document when it is no longer needed. Do not return it to the originator.

### Disclaimer

The views, opinions, and/or findings contained in this report are those of the author(s) and should not be construed as an official Department of the Army position, policy, or decision, unless so designated by other official documentation. Citation of trade names in this report does not constitute an official Department of the Army endorsement or approval of the use of such commercial items.

Reviewed:



---

JOHN V. BARSON  
LTC, MC, SFS  
Director, Biodynamics  
Research Division

Released for publication:



---

ROGER W. WILEY, O.D., Ph.D.  
Chairman, Scientific  
Review Committee



---

DAVID H. KARNEY  
Colonel, MC, SFS  
Commanding

## REPORT DOCUMENTATION PAGE

Form Approved  
OMB No. 0704-0188

1a. REPORT SECURITY CLASSIFICATION <b>UNCLASSIFIED</b>		1b. RESTRICTIVE MARKINGS	
2a. SECURITY CLASSIFICATION AUTHORITY		3. DISTRIBUTION / AVAILABILITY OF REPORT Approved for public release, distribution unlimited	
2b. DECLASSIFICATION / DOWNGRADING SCHEDULE		5. MONITORING ORGANIZATION REPORT NUMBER(S)	
4. PERFORMING ORGANIZATION REPORT NUMBER(S) USAARL Report No. 92-31		7a. NAME OF MONITORING ORGANIZATION U.S. Army Medical Research and Development Command	
6a. NAME OF PERFORMING ORGANIZATION U.S. Army Aeromedical Research Laboratory		7b. ADDRESS (City, State, and ZIP Code) Fort Detrick Frederick, MD 21702-5012	
6c. ADDRESS (City, State, and ZIP Code) P.O. Box 577 Fort Rucker, AL 36362-5292		9. PROCUREMENT INSTRUMENT IDENTIFICATION NUMBER Partial effort under contract No. DAMD 17-86-C-6215	
8a. NAME OF FUNDING / SPONSORING ORGANIZATION		8b. OFFICE SYMBOL (if applicable) SGRD-UAD-IE	
8c. ADDRESS (City, State, and ZIP Code)		10. SOURCE OF FUNDING NUMBERS	
		PROGRAM ELEMENT NO. 0603807A	PROJECT NO. 3M463807D83
		TASK NO. LC	WORK UNIT ACCESSION NO. 201
11. TITLE (Include Security Classification) Test and Evaluation Report of the Corometrics Medical Systems Neonatal Monitor, Model 506			
12. PERSONAL AUTHOR(S) James E. Bruckart, Martin Quattlebaum, Brendan Squire, and Bill Olding			
13a. TYPE OF REPORT Final		13b. TIME COVERED FROM _____ TO _____	
		14. DATE OF REPORT (Year, Month, Day) 1992 August	
		15. PAGE COUNT 61	
16. SUPPLEMENTARY NOTATION			
17. COSATI CODES			18. SUBJECT TERMS (Continue on reverse if necessary and identify by block number) Electromagnetic compatibility, test and evaluation, aeromedical equipment
FIELD	GROUP	SUB-GROUP	
01	03		
14	02		
19. ABSTRACT (Continue on reverse if necessary and identify by block number) The Corometrics Medical Systems Neonatal Monitor, Model 506, was tested for environmental and electromagnetic interference/compatibility in the UH-60A helicopter under the U.S. Army Program for Testing and Evaluation of Equipment for Aeromedical Operations. The tests were conducted using current military and industrial standards and procedures for environmental tests and electromagnetic interference/compatibility and human factors. The Corometrics Medical Systems Neonatal Monitor, Model 506, was found to be compatible with the U.S. Army MEDEVAC UH-60 Black Hawk.			
20. DISTRIBUTION / AVAILABILITY OF ABSTRACT <input checked="" type="checkbox"/> UNCLASSIFIED/UNLIMITED <input type="checkbox"/> SAME AS RPT. <input type="checkbox"/> DTIC USERS		21. ABSTRACT SECURITY CLASSIFICATION UNCLASSIFIED	
22a. NAME OF RESPONSIBLE INDIVIDUAL Chief, Scientific Information Center		22b. TELEPHONE (Include Area Code) (205) 255-6907	
		22c. OFFICE SYMBOL SGRD-UAX-SI	

## Table of contents

SECTION	PAGE
<b>1. <u>EXECUTIVE DIGEST</u></b>	
1.1 Test objectives . . . . .	1-1
1.2 Testing authority . . . . .	1-2
1.3 Scope . . . . .	1-2
1.4 Material description . . . . .	1-3
1.5 Summary . . . . .	1-3
1.6 Conclusion . . . . .	1-4
<b>2. <u>SUBTESTS</u></b>	
2.1 Initial inspection . . . . .	2-1
2.2 Battery life evaluation . . . . .	2-1
2.3 Electrical safety evaluation . . . . .	2-2
2.4 Human factors evaluation (laboratory) . . . . .	2-2
2.5 Altitude (low pressure) test . . . . .	2-3
2.6 Vibration test . . . . .	2-4
2.7 High temperature test . . . . .	2-5
2.8 Low temperature test . . . . .	2-7
2.9 Humidity test . . . . .	2-8
2.10 Electromagnetic characteristics test . . . . .	2-9
2.11 In-flight human factors evaluation . . . . .	2-11
2.12 In-flight EMI/EMC characteristics test . . . . .	2-12
<b>3. <u>SUPPORTING DOCUMENTATION</u></b>	
3.1 Detailed test information . . . . .	3-1
3.2 Test data . . . . .	3-2

3.3	Criteria, significant problems, and suggested improvements . . . . .	3-35
3.4	References . . . . .	3-38
3.5	Abbreviations . . . . .	3-39
3.6	List of manufacturers. . . . .	3-41
3.7	Distribution list . . . . .	3-42

## Section 1. Executive digest

The Army program for Test and Evaluation of Aeromedical Equipment uses existing military standards (MIL-STD) and collective professional expertise to test and evaluate selected medical equipment proposed for use aboard Army aircraft. Equipment meeting these standards ensures the safety of the crew, patients, and aircraft by eliminating risks due to: (1) Interference by the medical equipment with aircraft systems/subsystems operation, (2) the aircraft system's interference with the operation of the medical equipment, (3) the medical equipment's susceptibility to environmental exposure, or (4) physical and/or functional incompatibility while in use on board selected rotary-wing aircraft. This program tests both developmental and nondevelopmental (off the shelf) medical equipment destined for use aboard Army medical evacuation (MEDEVAC) aircraft.

### 1.1 TEST OBJECTIVES

1.1.1 To determine if the medical equipment is complete and operational per the manufacturer's operating instructions.

1.1.2 To ensure the electrical safety of the medical equipment.

1.1.3 To ensure the equipment will function as designed throughout the rated battery operation time.

1.1.4 To ensure the safety of the operator, the patient, and the aircrew.

1.1.5 To assess design considerations which could potentially contribute to an operator error.

1.1.6 To determine if the medical equipment can function as designed in a low pressure environment.

1.1.7 To determine the ability of the medical equipment to withstand the vibrational stresses expected in a rotary-wing flight environment without degradation or malfunction.

1.1.8 To determine the ability of the medical equipment to be stored and operated in a high temperature environment.

1.1.9 To determine the ability of the medical equipment to be stored and operated in a low temperature environment.

1.1.10 To determine the ability of the medical equipment to operate satisfactorily for short periods during exposure to high humidity conditions.

1.1.11 To assess the levels of electromagnetic emissions produced by the medical equipment within selected frequency ranges.

1.1.12 To assess the minimum electromagnetic susceptibility levels of the medical equipment within selected frequency ranges.

1.1.13 To assess the physical and/or functional compatibility of the medical equipment while in use on board the aircraft.

1.1.14 To assess the electromagnetic interference (EMI) and electromagnetic compatibility (EMC) characteristics of the medical equipment with the host aircraft and its installed systems.

## 1.2 TESTING AUTHORITY

Research and Technology Work Unit Summary, dated 5 October 1989. Project number 3M463807D836, titled, Army Program for Testing and Evaluation of Equipment for Aeromedical Operations.

## 1.3 SCOPE

1.3.1 This test was conducted at the United States Army Aeromedical Research Laboratory (USAARL), Cairns Army Airfield (CAAF), and designated test flight areas in and around Fort Rucker, Alabama.

1.3.2 The USAARL UH-60A aircraft, serial number 88-26069, with subsystems delineated in paragraph 3.2.2, was configured with the Corometrics Medical Systems Neonatal Monitor\*, model 506 and used as the test aircraft for the in-flight evaluation. The in-flight evaluation required 3.3 flight hours.

1.3.3 Laboratory testing was accomplished at USAARL using government furnished equipment (GFE) by Universal Energy Systems, Inc. (UES), under contract No. DAMD 17-86-C-6215.

1.3.4 Prior to flight testing, the following tests were accomplished: Acceptance inspection, equipment training, electromagnetic compatibility, human factors and safety, environmental compatibility, and in-flight compatibility.

1.3.5 An airworthiness release (AWR) dated 24 Feb 1992 was received from the U.S. Army Aviation Systems Command (AVSCOM) prior to the in-flight testing of the Corometrics Neonatal Monitor.

---

\* See list of manufacturers

## 1.4 MATERIAL DESCRIPTION

The Corometrics Neonatal Monitor, Model 506, is designed to monitor newborn infants during intensive care and transport. The monitor displays the infant's heart rate, respiration rate, systolic and diastolic blood pressure, and two temperatures. The electrocardiogram (ECG), respiration, and blood pressure are shown on a cathode ray tube (CRT). Sliding levers on the front panel allow selection of upper and lower limits for heart rate and delay time between breaths (apnea) or maximum number of breaths (tachypnea). A switch on the front of the unit allows selection of TEMP 1, TEMP 2, or MEAN PRES (mean pressure) display. Adjusting knobs allow variation in the size of the ECG trace and respiration trace. A "DISPLAY" button applies power to the monitor and "ALARMS OFF" disables the audible alarm. "FREEZE" stops the trace on the CRT.

The Corometrics Neonatal Monitor uses a microprocessor to allow continuous monitoring of ECG, heart rate, blood pressure, and respiration rate. ECG, heart rate, and respiration rate are derived from the ECG signal. Blood pressure is measured from a pressure transducer. Temperature is measured with a thermistor. The unit operates from ac power or an internal rechargeable NiCad battery.

## 1.5 SUMMARY

### 1.5.1 Laboratory testing

1.5.1.1 Battery Life Evaluation: During three battery life cycle tests, the unit averaged 5 hours and 53 minutes of operation at room temperature while monitoring a simulated ECG signal. This exceeds the manufacturer specification of 4 hours operation on a fully charged battery.

1.5.1.2 Electrical Safety Evaluation: All measurements were within acceptable limits. No unsafe qualities were found in the Corometrics Model 506. The limits for currents and resistances were in accordance with (IAW) the limits specified in TB-38-750-2, April 1987 and National Fire Prevention Association (NFPA) standards.

1.5.1.3 Human Factors Evaluation: The Corometrics Model 506 was found to be satisfactory in all except four categories of the evaluation. This unit uses the "DISPLAY" switch to control power instead of "ON/OFF" or "POWER". The CRT display does not have an illumination control. Switch and knob identification/position labels are not illuminated. The unit does not have an audible alarm for a low battery condition.

1.5.1.4 Environmental Tests: The Corometrics Model 506 can be expected to perform in a variety of environmental conditions.

Its performance was found to be satisfactory in all stages of the environmental testing. The requirements for environmental tests are established in MIL-STD-810D, Methods 500.2 (altitude), 514.3 (vibration), 501.2 (high temperature), 502.2 (low temperature), and 507.2 (humidity).

1.5.1.5 Radiated Emissions Tests (RE02): The Corometrics Model 506 may be unsatisfactory for use in certain EMI sensitive environments. Broadband (BB) and narrowband radiated emissions were detected in the test frequency ranges. Some emissions exceeded the test limits. Emission limits are set forth in MIL-STD-461A, Notice 4.

1.5.1.6 Radiated Susceptibility Test (RS03): The Corometrics Model 506 was susceptible to radio frequency interference in the testing range and magnitude. Failure was noted as erratic wave forms, inaccurate numeric displays, and erroneous alarms.

1.5.1.7 Conducted Emissions Test (CE01, CE02, and CE04): No signal failures were detected from the Corometrics Neonatal Monitor during this test.

1.5.1.8 Conducted Susceptibility Test (CS02 and CS06): No susceptibility to the test power line spikes was noted in the Corometrics Monitor.

#### 1.5.2 In-flight testing

1.5.2.1 During the in-flight human factors evaluation, the Corometrics Model 506 was found to be satisfactory in all categories of the evaluation criteria. Vibration of the ECG leads caused some distortion of the ECG waveform on the CRT. This distortion was not present if the leads were not allowed to vibrate with the aircraft structure.

1.5.2.2 The aircraft and its subsystems were not adversely affected by the operation of the Corometrics Model 506 in any of the prescribed flight test modes.

1.5.2.3 The Corometrics Model 506 was not affected by the aircraft and its subsystems during the in-flight testing.

#### 1.6 CONCLUSION

Based on the results of laboratory and in-flight testing, the Corometrics Medical Systems Neonatal Monitor, Model 506 was found to be compatible with U.S. Army MEDEVAC UH-60A Black Hawk with the subsystems listed in paragraph 3.2.2.

## Section 2. Subtests

### 2.1 INITIAL INSPECTION

#### 2.1.1 Objective

To determine if the Corometrics Model 506 is complete and operational for testing per the manufacturer's operating instructions.

#### 2.1.2 Criteria

2.1.2.1 The physical inventory is conducted solely for investigation and documentation.

2.1.2.2 The Corometrics Model 506 will display consistent and accurate performance as an acceptable performance test.

#### 2.1.3 Test procedure

2.1.3.1 A complete physical inventory of the Corometrics Model 506 was completed per the manufacturer's equipment list.

2.1.3.2 An operational validation test of the Corometrics Model 506 was conducted per the manufacturer's operating instructions by USAARL's medical maintenance personnel.

#### 2.1.4 Test findings

2.1.4.1 The Corometrics Model 506 was inventoried and found to be complete.

2.1.4.2 The Corometrics Model 506 operated as prescribed in the manufacturer's operating manual. Criteria met.

### 2.2 BATTERY LIFE EVALUATION (Laboratory)

#### 2.2.1 Objective

To ensure the equipment will function as designed throughout the rated battery operation time.

#### 2.2.2 Criterion

Verify manufacturer's specified full power internal battery life expectancy of 4 hours operation.

#### 2.2.3 Test procedure

2.2.3.1 Charging and operation cycles were conducted in ambient room conditions.

#### 2.2.4 Test findings

The battery provided 5 hours and 53 minutes service while monitoring a simulated ECG signal. This exceeds the manufacturers specification for battery life. Criterion met.

### 2.3 ELECTRICAL SAFETY EVALUATION

#### 2.3.1 Objective

To ensure the electrical safety, by evaluation of case-to-ground resistance and case-to-ground current leakage, of the Corometrics Model 506.

#### 2.3.2 Criterion

The Corometrics Model 506 shall meet the standards established in TB-38-750-2 and NFPA 99 for electrical safety of medical equipment.

#### 2.3.3 Test procedure

Performance in the electrical safety evaluation were made, with a Neurodyne-Dempsey model 431F electrical safety analyzer\*, IAW the procedures described in Technical Bulletin (TB) Number 38-750-2. Case-to-ground resistance and various case-to-ground leakage currents were measured. Leakage currents were measured using a 10 by 20 centimeter (cm) aluminum foil sheet taped flush to the equipment case. Checks were made for safety concerns such as case integrity, breaks in power cord insulation, and connectors.

#### 2.3.4 Test findings

Grounding conductor resistance was 63 milliohms, maximum case leakage current was 18.8 microamperes, and maximum lead leakage current was 3.0 microamperes. These measurements are below the limits specified in TB-38-750-2 and NFPA 99. Criterion met.

### 2.4 HUMAN FACTORS EVALUATION (Laboratory)

#### 2.4.1 Objectives

2.4.1.1 To assure the safety of the operator, the potential patient, and the aircrew.

2.4.1.2 To assess the design considerations which could potentially contribute to an operator error.

#### 2.4.2 Criterion

The Corometrics Model 506 must be rated satisfactory in all major categories of the evaluation. These include visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.

#### 2.4.3 Test procedure

2.4.3.1 The evaluation was conducted in a laboratory under fluorescent lighting and ambient room conditions.

2.4.3.2 The Corometrics Model 506 was operated according to prescribed instructions through its full range of functions.

#### 2.4.4 Test finding

The first Corometrics Model 506 was found to be satisfactory in all of the evaluation criteria except Controls and Labels and Coding. The unit's power switch is labeled "DISPLAY" instead on "ON/OFF" or "POWER." The CRT display does not have an illumination control. Switch and knob position labels are not illuminated. The unit does not have an audible alarm to indicate a low battery condition. Criterion partially met.

### 2.5 ALTITUDE (LOW PRESSURE) TEST [IAW MIL-STD-810D, METHOD 500.2]

#### 2.5.1 Objective

To determine if the Corometrics Model 506 can function as designed in a low pressure environment.

#### 2.5.2 Criterion

The Corometrics Model 506 will perform as designed while exposed to an altitude equivalency of 15,000 feet above sea level.

#### 2.5.3 Test procedure

2.5.3.1 A pretest performance check was conducted to ensure proper operation of the Corometrics Model 506.

2.5.3.2 The altitude test was performed in a Tenney Engineering model 64S altitude chamber\*. This test is based on MIL-STD-810D, Method 500.2. The Corometrics Model 506 was operated on the

floor of the chamber. Chamber pressure was decreased to 420 mmHg (15,000 ft equivalent altitude) over a 15-minute period, held constant for 60 minutes, then raised, at 1500 fpm, to ambient conditions (760 mmHg) over a 10-minute period. There are no provisions for the control of temperature or humidity inside this chamber.

2.5.3.3 A posttest performance check was conducted to ensure proper operation of the Corometrics Model 506 after the exposure to low pressure.

#### 2.5.4 Test findings

2.5.4.1 The pretest performance check met criterion 2.1.2.2.

2.5.4.2 No failures in the performance of the Corometrics Model 506 were noted before, during, or after the altitude test. Criterion met.

2.5.4.3 The posttest performance check met criterion 2.1.2.2.

### 2.6 VIBRATION TEST [IAW MIL-STD-810D, METHOD 514.3]

#### 2.6.1 Objective

To determine the ability of the Corometrics Model 506 to withstand the vibrational stresses expected in a rotary-wing environment without degradation or malfunction.

#### 2.6.2 Criterion

The Corometrics Model 506 will remain operational and be able to display consistent and accurate performance while exposed to vibrational stresses.

#### 2.6.3 Test procedure

2.6.3.1 A pretest performance check was conducted to ensure proper operation of the Corometrics Model 506.

2.6.3.2 The vibration test was performed using an Unholtz-Dickey model TA115-40/CSTA vibration test system\*. It is a single-axis system with an electromagnetic driver unit. The test consisted of sinusoidal vibrations superimposed on random vibrations over a frequency range of 500 Hz, as shown below. These vibrations are derived from performance taken on the floor under the copilot's seat in a UH-1 helicopter traveling at 120 knots. The reference spectrum breakpoints are from MIL-STD-810D, Method 514.3; reference spectrum levels are based on field performance with a conservatism factor of 1.5. Independent tests were conducted in the X, Y, and Z axes.

### Z-axis

duration: 60 minutes  
broadband intensity: 0.4506  $G_{rms}$   
random vibration: initial slope : 99.00 dB/oct  
5 Hz level: 0.00006210  $G_{sqr/Hz}$   
100 Hz level: 0.0006210  $G_{sqr/Hz}$   
300 Hz level: 0.0006210  $G_{sqr/Hz}$   
500 Hz level: 0.00006210  $G_{sqr/Hz}$   
final slope: -99.00 dB/oct  
sinusoidal vibration: .5450  $G_{pk}$  at 11.25 Hz  
.1690  $G_{pk}$  at 22.50 Hz  
.1200  $G_{pk}$  at 33.75 Hz  
.0310  $G_{pk}$  at 45.00 Hz  
.0530  $G_{pk}$  at 56.25 Hz

### X and Y axes

duration: 60 minutes each  
broadband intensity: 0.3099  $G_{rms}$   
random vibration: initial slope: 99.00 dB/oct  
5 Hz level: 0.00002920  $G_{sqr/Hz}$   
100 Hz level: 0.0002920  $G_{sqr/Hz}$   
300 Hz level: 0.0002920  $G_{sqr/Hz}$   
500 Hz level: 0.00002920  $G_{sqr/Hz}$   
final slope: -99.00 dB/oct  
sinusoidal vibration: .3200  $G_{pk}$  at 11.25 Hz  
.0670  $G_{pk}$  at 22.50 Hz  
.0950  $G_{pk}$  at 33.75 Hz  
.0350  $G_{pk}$  at 45.00 Hz  
.0770  $G_{pk}$  at 56.25 Hz

The Corometrics Model 506 was strapped to the vibration table fixture, and its performance was evaluated before, during, and after exposure to vibration.

2.6.3.3 A posttest performance check was conducted to ensure proper operation of the Corometrics Model 506.

### 2.6.4 Test findings

2.6.4.1 The pretest performance check met criterion 2.1.2.2.

2.6.4.2 No failures in the performance of the Corometrics Model 506 occurred before, during, or after exposure to vibration. Criterion met.

2.6.4.3 The posttest performance check met criterion 2.1.2.2.

## 2.7 HIGH TEMPERATURE TEST [IAW MIL-STD-810D, METHOD 501.2]

### 2.7.1 Objective

To determine the ability of the Corometrics Model 506 to be stored and operated in a high temperature environment.

### 2.7.2 Criteria

2.7.2.1 The Corometrics Model 506 will demonstrate consistent and accurate operation during the high temperature operation check.

2.7.2.2 The Corometrics Model 506 will demonstrate consistent and accurate operation after the high temperature storage cycle.

### 2.7.3 Test procedure

2.7.3.1 A pretest performance check was conducted to ensure proper operation of the Corometrics Model 506.

2.7.3.2 The high temperature test was conducted in a Tenney Engineering model ZWUL-10107D walk-in controlled environment chamber\*. This test is based on MIL-STD-810D, Method 501.2. For the high temperature operation test, the Corometrics Model 506 was turned on and placed on the floor of the environmental chamber. The chamber temperature was raised to 49°C and the humidity was stabilized at a maximum of 20 percent relative humidity (RH) within 15 minutes. The environmental control system is capable of regulating temperature within  $\pm 2^\circ\text{C}$  and humidity within  $\pm 5$  percent RH. Temperature and humidity were held constant for 2 hours. At 30-minute intervals, the chamber door was opened briefly to minimize the change in chamber conditions during performance checks. After the operational test, the Corometrics Model 506 was allowed to return to ambient conditions over a 30-minute period.

2.7.3.3 A posttest performance check was conducted to ensure proper operation of the Corometrics Model 506.

2.7.3.4 The Corometrics Model 506 was stored (not operated) at temperatures of 63°C for 1 hour, 71°C for 4 hours, then again at 63°C for 1 hour. The chamber and Corometrics Model 506 then were returned to ambient conditions over a 30-minute period.

2.7.3.5 A poststorage performance check was conducted to ensure proper performance of the Corometrics Model 506.

### 2.7.4 Test findings

2.7.4.1 The pretest performance check met criterion 2.1.2.2.

2.7.4.2 No operational failures occurred during the high temperature test. Criterion met.

2.7.4.3 The posttest performance check met criterion 2.1.2.2.

2.7.4.4 The Corometrics Model 506 functioned properly after the high temperature storage test. Criterion met.

## 2.8 LOW TEMPERATURE TEST [IAW MIL-STD-810D, METHOD 502.2]

### 2.8.1 Objective

To determine the ability of the Corometrics Model 506 to be stored and operated in a low temperature environment.

### 2.8.2 Criteria

2.8.2.1 The Corometrics Model 506 will demonstrate consistent and accurate operation during the low temperature operation check.

2.8.2.2 The Corometrics Model 506 will demonstrate consistent and accurate operation after the low temperature storage cycle.

### 2.8.3 Test procedure

2.8.3.1 A pretest performance check was conducted to ensure proper operation of the Corometrics Model 506.

2.8.3.2 The Corometrics Model 506 was placed on the floor of the environmental chamber and the temperature was lowered to 0°C within 25 minutes. The environmental control system is capable of regulating temperature within 2°C. Humidity cannot be controlled in the chamber at freezing temperatures. The temperature was held constant for 2 hours. The chamber door was opened briefly every 30 minutes to minimize the change in chamber conditions, and a performance check was conducted. The chamber temperature then was raised to ambient temperature within a 30-minute period.

2.8.3.3 A posttest performance check was conducted to ensure proper operation of the Corometrics Model 506.

2.8.3.4 The Corometrics Model 506 was "stored" in a nonoperational mode. The Corometrics Model 506 was placed on the floor of the environmental test chamber and the temperature was lowered to -46°C for 6 hours. The chamber then was raised to ambient temperature over a 30-minute period.

2.8.3.5 A poststorage performance check was conducted to ensure proper operation of the Corometrics Model 506.

### 2.8.4 Test findings

2.8.4.1 The pretest performance check met criterion 2.1.2.2.

2.8.4.2 No operational failures occurred during the low temperature test. Criterion met.

2.8.4.3 The posttest performance check met criterion 2.1.2.2.

2.8.4.4 The Corometrics Model 506 functioned properly after the low temperature storage test. Criterion met.

## 2.9 HUMIDITY TEST [IAW MIL-STD-810D, METHOD 507.2]

### 2.9.1 Objective

To determine the ability of the Corometrics Model 506 to operate satisfactorily for short periods of time during exposure to highly humid conditions.

### 2.9.2 Criterion

The Corometrics Model 506 will demonstrate consistent and accurate operation while exposed to a high humidity environment.

### 2.9.3 Test procedure

2.9.3.1 A pretest performance check was conducted to ensure the proper operation of the Corometrics Model 506.

2.9.3.2 The humidity test was conducted in a Tenney Engineering model ZWUL-10107D walk-in controlled environment chamber\*. This test is based on MIL-STD-810D, Method 507.2. For the humidity test, the Corometrics Model 506 was placed in operation on the floor of the environmental chamber. The chamber temperature was raised to a temperature of 30°C and a relative humidity of 95 percent within 25 minutes. Temperature and relative humidity were maintained for 4 hours. The environmental control system is capable of regulating temperature within  $\pm 2^\circ\text{C}$  and humidity within  $\pm 5$  percent RH. At 45-minute intervals the performance of the Corometrics Model 506 was checked. The chamber door was opened briefly to minimize the change in chamber conditions. The chamber and the Corometrics Model 506 were returned to ambient conditions before the posttest performance validation check was conducted.

2.9.3.3 A posttest performance check was conducted to ensure the proper operation of the Corometrics Model 506.

### 2.9.4 Test findings

2.9.4.1 The pretest performance check met criterion 2.1.2.2.

2.9.4.2 No failures were noted in the Corometrics Model 506 performance checks conducted during the exposure to the high humidity environment. Criterion met.

2.9.4.3 The posttest performance check met criterion 2.1.2.2.

2.10 ELECTROMAGNETIC CHARACTERISTICS TEST [IAW MIL-STD-461A, Notice 4, AND MIL-STD-462, Notice 3]

#### 2.10.1 Objectives

2.10.1.1 To assess the maximum levels of radiated electromagnetic emissions produced by the Corometrics Model 506 in the 14 kHz to 12.4 GHz frequency range.

2.10.1.2 To assess the tolerances of radiated electromagnetic susceptibility of the Corometrics Model 506 within the 10 kHz to 10 GHz electric field.

2.10.1.3 To assess the maximum levels of conducted electromagnetic emissions produced by the Corometrics Model 506 in the 10 kHz to 50 MHz frequency ranges.

2.10.1.4 To assess the tolerances of conducted electromagnetic susceptibility of the Corometrics Model 506 within the range of 50 kHz to 400 MHz and power spikes.

#### 2.10.2 Criteria

2.10.2.1 The Corometrics Model 506 will not produce emissions in excess of the limits set forth in MIL-STD-461A, Notice 4, paragraph 6.13.

2.10.2.2 The Corometrics Model 506 will not malfunction when it is subjected to radiated emissions as specified in MIL-STD-461A, Notice 4, paragraph 6.20.

2.10.2.3 The Corometrics Model 506 will not conduct emissions in excess of the limits set forth in MIL-STD-461A, Notice 4, paragraphs 6.1 and 6.2.

2.10.2.4 The Corometrics Model 506 will not malfunction when it is subjected to conducted emissions as specified in MIL-STD-461A, Notice 4, paragraphs 6.7 and 6.10.

#### 2.10.3 Test procedure

2.10.3.1 The radiated emissions test was performed according to MIL-STD-462, Notice 3, Method RE02. The Corometrics Model 506 was positioned on a wooden test stand inside the EMI chamber, 1 meter away from the receiving antennas. The antennas were

mounted for both vertical and horizontal polarities and connected to EMI receivers. While the Corometrics Model 506 was operating, the frequency spectrum (14 kHz to 12.4 GHz) was scanned for emissions. The Corometrics Model 506 was operated with ac and battery power.

2.10.3.2 The radiated susceptibility test was performed according to MIL-STD-462, Notice 3, Method RS03. The Corometrics Model 506 was positioned on a wooden test stand inside the EMI chamber 1 meter away from the transmitting antennas. The antennas were mounted for both vertical and horizontal polarities and connected to radio frequency (RF) transmitters. While the Corometrics Model 506 was operating, it was monitored for faulty operation during exposures to fields of 1 V/m from 10 kHz to 2 MHz, and 5 V/m from 2 to 30 MHz, 10 V/m from 30 MHz to 2 GHz, and 5 V/m from 2 to 10 GHz. The Corometrics Model 506 was operated with ac and battery power.

2.10.3.3 The conducted emissions tests were performed according to MIL-STD-462, Notice 3, Methods CE02 and CE04. The Corometrics Model 506 was placed on a grounded, copper-covered workbench. The top of the workbench was 1 meter from floor level, 1.37 meters long and 0.81 meters wide. Power was supplied via a pair of line impedance stabilization networks (LISN) and a test jig. The test jig is a wooden tray with two power receptacles and two slots to hold current probes in place around power supply conductors. While the Corometrics Model 506 was operating, the frequency range (10 kHz to 50 MHz) was scanned for emissions conducted in the power cable from the Corometrics Model 506.

2.10.3.4 The conducted susceptibility spike test was performed on a chemical resistant counter top according to MIL-STD-462, Notice 3, Method CS06. Power was supplied via a customized metal connection box. The connection box has two power receptacles and four banana jacks on its front panel. Connections to the individual power lines were made in series through the banana jacks. Transient spikes of 100 volts, 10 microseconds were generated with a Solar Electronics model 8282-1 transient pulse generator\* and induced onto the power leads at the connection box banana jacks. The spikes were monitored with a Tektronix 2235 oscilloscope\* connected to a power receptacle on the connection box. The Corometrics Model 506 was plugged into the other receptacle on the connection box and placed in operation. It was observed visually for correct operation while it was subjected to the power line spikes.

2.10.3.5 The conducted susceptibility test was performed according to MIL-STD-462, Notice 3, Method CS02. The Corometrics Model 506 was placed on a grounded, copper-covered workbench. Radio frequency interference was induced on the power leads and measured at the Corometrics Model 506 power cable. The frequency of

the interference was incremented over the 50 kHz to 400 MHz range while the Corometrics Model 506 was operated. It was observed visually for proper operation while it was subjected to the radio interference on the power leads. Each frequency was held for 15 seconds.

#### 2.10.4 Test findings

2.10.4.1 During the radiated emissions test, emissions which exceeded specification limits of MIL-STD-461A, Notice 4, were detected. These included:

	<u>Frequency range</u>	<u>Emission exceeding standard</u>
ac:	0.066 - 107.2 MHz	0.2 - 25.5 dB (NB)
	0.014 - 0.15 MHz	3.2 - 10.4 dB (BB)
battery:		
	0.015 - 107.275 MHz	0.4 - 36.5 dB (NB)
	0.014 - 0.55 MHz	4.3 - 25.1 dB (BB)
	1.025 MHz	1.3 dB (BB)
	1.672 MHz	8.5 dB (BB)

Criterion partially met.

2.10.4.2 The Corometrics Model 506 was susceptible to radio frequency interference in the testing range and magnitude.

	<u>Frequency</u>	<u>Field strength</u>
ac:	6.7 - 7.7 MHz	0.53 - 1.67 V/m
	14.8 - 22.4 MHz	0.28 - 3.34 V/m
	30.0 - 176.2 MHz	0.63 - 7.94 V/m
battery:		
	20.4 - 24.8 MHz	0.35 - 4.21 V/m
	30.0 - 183 MHz	0.53 - 8.91 V/m

Criterion partially met.

2.10.4.3 No signal failures were detected from the Corometrics Model 506 during the conducted emissions test. Criterion met.

2.10.4.4 The Corometrics Model 506 was not susceptible to radio frequency interference (RFI) or test spikes during the conducted susceptibility tests. Criterion met.

#### 2.11 IN-FLIGHT HUMAN FACTORS EVALUATION

##### 2.11.1 Objective

To assess the physical and/or functional compatibility of the Corometrics Model 506 while in use onboard the aircraft.

## 2.11.2 Criterion

The flight surgeon will be able to operate the Corometrics Model 506 without physical or functional restrictions aboard the aircraft. Major areas of concern include: Proper operation, visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.

## 2.11.3 Test procedure

2.11.3.1 A human factors evaluation was performed IAW MIL-STD-1472D, AAMI Human factors engineering guidelines, and UL-544 to ensure the compatibility of the Corometrics Model 506 and the in-flight environment. The flight surgeon conducted the test wearing a flight suit, flight gloves, and an SPH-4B flight helmet. An evaluation of the compatibility with the nuclear, biological, and chemical (NBC) protective equipment was not conducted. Due to restrictions of the AWR, testing was conducted during daylight hours only.

2.11.3.2 The Corometrics Model 506 was placed on the floor of the aircraft and secured with cargo straps. The Corometrics Model 506 was tested using ac and battery power in all flight scenarios required by the In-flight Test Operations Procedures (ITOP).

## 2.11.4 Test findings

During the in-flight human factors evaluation, the Corometrics Model 506 was found to be satisfactory in all categories of the evaluation criteria except as noted in the laboratory human factors evaluation. When the patient leads vibrated in response to aircraft vibration, some interference was seen in the baseline ECG signal. The high frequency "noise" in the ECG signal disappears if the leads are secured by taping to a fixed object. Criterion partially met.

## 2.12 IN-FLIGHT EMI/EMC CHARACTERISTICS

### 2.12.1 Objective

To assess the EMI/EMC characteristics of the Corometrics Model 506 with the host aircraft and its installed systems.

### 2.12.2 Criteria

2.12.2.1 The Corometrics Model 506 will not radiate EMI to disrupt or interfere with other equipment or systems aboard the aircraft.

2.12.2.2 The aircraft will not radiate EMI to disrupt or interfere with the Corometrics Model 506's operation.

2.12.3 Test procedure

A qualitative EMI/EMC assessment was performed with both the Corometrics Model 506 and the aircraft operating as source and victim. The Corometrics Model 506 and applicable aircraft instruments and systems were monitored for unusual operation, readings, surges, or power anomalies for each checklist item (see section 3).

2.12.4 Test findings

2.12.4.1 There were no adverse instances of EMI/EMC noted with the Corometrics Model 506 acting as either the source or victim. Criterion met.

2.12.4.2 There were no adverse instances of EMI/EMC noted with the aircraft acting as either the source or victim. Criterion met.

## Section 3. Supporting documentation

### 3.1 DETAILED TEST INFORMATION

#### 3.1.1 General information

3.1.1.1 Corometrics Model 506 testing is not considered a major action significantly affecting the quality of the human environment and, therefore, qualifies for categorical exclusion A-28, appendix A, AR 200-1.

3.1.1.2 A safety pilot will be designated for each flight. Flight operations will be conducted IAW the aircraft operator's manual, appropriate aircrew training manuals, and test item technical data.

#### 3.1.2 Material description

3.1.2.1 The Corometrics Neonatal Monitor, Model 506, is designed to monitor newborn infants during intensive care and transport. The monitor displays the infant's heart rate, respiration rate, systolic and diastolic blood pressure, and two temperatures. The electrocardiogram (ECG), respiration, and blood pressure are shown on a cathode ray tube (CRT). Sliding levers on the front panel allow selection of upper and lower limits for heart rate and delay time between breaths (apnea) or maximum number of breaths (tachypnea). A switch on the front of the unit allows selection of TEMP 1, TEMP 2, or MEAN PRES (mean pressure) display. Adjusting knobs allow variation in the size of the ECG trace and respiration trace. A "DISPLAY" button applies power to the monitor and "ALARMS OFF" disables the audible alarm. "FREEZE" stops the trace on the CRT.

The Corometrics Neonatal Monitor uses a microprocessor to allow continuous monitoring of ECG, heart rate, blood pressure, and respiration rate. ECG, heart rate, and respiration rate are derived from the ECG signal. Blood pressure is measured from a pressure transducer. Temperature is measured with a thermistor. The unit operates from ac power or an internal rechargeable NiCad battery.


3.1.2.2 Dimensions: 15.2 x 23.9 x 36.6 cm (4 x 9.4 x 14.4 in).

3.1.2.3 Weight: 7.9 kg (17.5 lbs)

3.1.2.4 Power requirements: 117 Vac, 50/60 Hz, 20 watts. An internal, rechargeable NiCad battery provides a rated 4 hours of operation on a fully charged battery. Charge time for a completely discharged battery to full charge is 14 to 16 hours.

3.2 TEST DATA

3.2.1 Photographic description

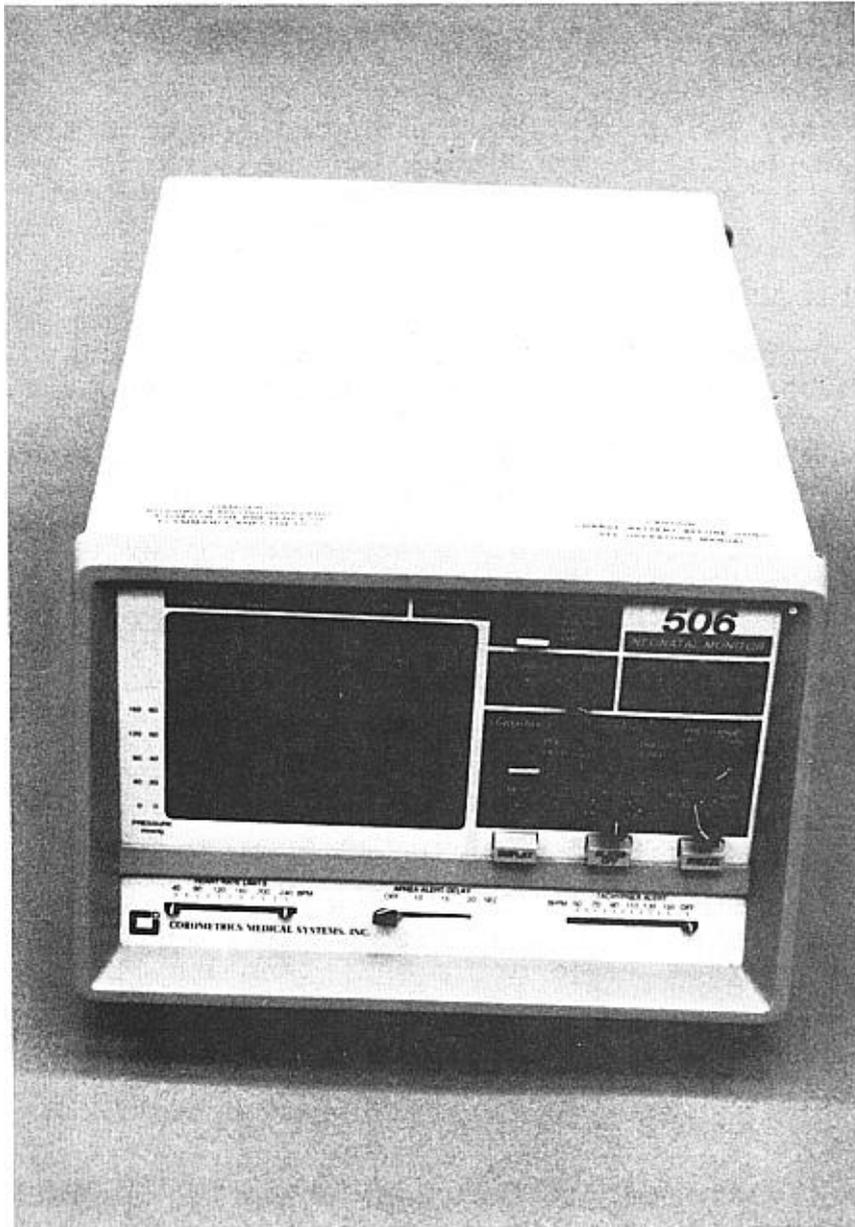

### 3.2.2 Aircraft equipment list

Item No.	Nomenclature
1	Receiver radio -- R-1496A/ARN-89 (automatic direction finder)
2	Displacement gyro -- CN-1314/A
3	Gyro directional -- CN-998/ASN-43
4	Signal data converter -- CV-3338/ASN-128
5	Receiver -- R-2139/ARN-123 (VOR/LOC/MB/GS)
6	Command instrument system processor -- 70600-01038-101
7	SAS amplifier -- 70901-02908-104 (flight control stability augmentation system)
8	Rate gyro -- TRU-2A/A
9	Amplifier, impedance -- AM-4859A/ARN-89
10	Cargo hook -- FE-7590-145
11	Receiver, radar -- RT-1193/ASN-128 (doppler navigation receiver)
12	Barometric altimeter -- AAU-31/A-1
13	Barometric altimeter -- AAU-32A
14	Receiver/transmitter -- RT-1300/ARC-186 (VHF-AM and/or FM radio)
15	UHF-AM radio set -- RT-1518/ARC-164
16	Interphone control -- C6533/ARC (aircraft intercom control)
17	Receiver/transmitter -- RT-1115D/APN-209 (radar altimeter)
18	Indicator altimeter -- ID-1917C/APN-209 (radar altimeter)
19	Control radio set -- C-7392A/ARN-89 (automatic direction finder)
20	Comparator signal data -- CM-482/ARC-186 (comparator for ARC-186)
21	Receiver/transmitter -- RT-1296A/APX-100 (transponder with IFF)
22	Computer display unit -- CP-1252/ASN-128 (doppler navigation system)
23	Compass set controller -- C-8021E/ASN75
24	Magnetic compass - standby -- MS-17983-4

### 3.2.3 In-flight test data card

#### DATA CARD FORMAT

#### GUIDELINE FOR DATA COLLECTION

#### IN-FLIGHT SUITABILITY TEST OF MEDICAL ITEMS

1. Installation/removal.	Suitable		Comments
	Yes	No	
a. Weight and balance (DD Form 365-4, Clearance Form F).	X		
b. Space/area allocation.			
(1) Operational requirements.	X		
(2) Storage requirements.	X		
c. Interface connections (safe, positive, secure).	X		
d. Installation/removal (expedient/easily achieved).	X		
e. Mounting/final config- uration (functional/stable).	X		
2. Operations and performance.	Suitable		Comments
	Yes	No	
a. Manufacturer's operating instruction.	X		
b. Medical item operation before aircraft run-up.	X		
c. System interface during aircraft engine run-up and medical item operation (EMI switchology checklist).	X		
(1) Aircraft voltage output.	X		

	Suitable		Comments
	Yes	No	
(2) Flight control function (UH-60).	X		
(3) Stabilator function (UH-60).	X		
(4) Radio communication vs. medical item operation.			
(a) FM	X		
(b) UHF	X		
(c) VHF	X		
(5) Navigation equipment vs. medical item operation.			
(a) Transponder	X		
(b) ADF	X		
(c) VOR	X		
(d) Doppler	X		
(6) Radar altimeter operation vs. medical item operation.	X		
d. System interface during aircraft hover and medical item operation (EMI switchology checklist).			
(1) Voltage output.		NA	
(2) Radio communication vs. medical item operation.			
(a) FM	X		
(b) UHF	X		
(c) VHF	X		

(3) Navigation equipment operation vs. medical item operation.	Suitable		Comments
	Yes	No	
(a) Transponder	X		
(b) ADF	X		
(c) VOR	X		
(d) Doppler	X		

e. Flight mission profile vs. medical item operation (EMI switchology checklist).

(1) Straight and level (1000 ft MSL for 20 minutes).

(a) Compatibility of flight mode and medical item operation. X

(b) Radio communication vs. medical item operation.

a. FM X

b. UHF X

c. VHF X

(2) NOE (20 minutes). compatibility of flight mode and medical item operation. X

(3) FM homing (10 minutes). X

(4) Doppler navigation vs. medical item operation.

(a) Initialize function. X

(b) Fix function. X

(c) Update function. X

	Suitable		Comments
	Yes	No	
(5) VOR navigation 7000 ft MSL for 20 minutes) vs. medical item operation.	X		
(6) ILS approach vs. medical item operation.	X		
f. Medical item operation after engine shutdown (external power source).	X		
g. Restrictions to the medical item's use (i.e., electrical connectors).	X		
h. Deviations from the labor- atory test results.			
(1) Electrical/ electronic.		None	
(2) Mechanical environment.		None	
(3) Human factors (user interface, controls, markings, lighting, egress).		None	
(4) Safety.		None	
3. Deviations from the in-flight test protocol.			
a. The VOR navigation portion of the in-flight test con- ducted at 2000 feet MSL due to air traffic control clearance.			

3.2.4 EMI switchology checklist

EMI SWITCHOLOGY CHECKLIST UH-60 AIRCRAFT  
IN-FLIGHT SUITABILITY OF MEDICAL ITEMS

ENG INSTRUMENTS/CDU	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
Fuel quantity	X			
Fuel indicator test	X			
XMSN oil temperature	X			
XMSN oil pressure	X			
#1 engine oil temperature	X			
#2 engine oil temperature	X			
#1 engine oil pressure	X			
#2 engine oil pressure	X			
#1 TGT	X			
#2 TGT	X			
#1 Ng speed	X			
#2 Ng speed	X			
CDU digits on/off	X			
CDU instruments dim	X			

ENG INSTRUMENTS/PLT PDU	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
#1 engine RPM	X			
#2 engine RPM	X			
Rotor RPM	X			
#1 torque	X			
#2 torque	X			

ENG INSTRUMENTS/COPLT PDU	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
#1 engine RPM	X			
#2 engine RPM	X			
Rotor RPM	X			
#1 torque	X			
#2 torque	X			

ENG CONTROLS	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
#1 overspeed	X			
#2 overspeed	X			
RPM switch	X			
#1 engine anti-ice	X			
#2 engine anti-ice	X			
#1 inlet anti-ice	X			
#2 inlet anti-ice	X			

RADIO EQUIPMENT	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
ICS, C-6533 ARC	X			
VHF-FM, ARC-186/115	X			
VHF-AM, ARC-186/115	X			
UHF-AM, ARC-164(V)	X			
Crypto, KY-28	Not installed			
Radio retransmissions PLN	Not installed			
Transponder, APX-100(V)	X			
KIT-1A/TSEC IFF computer	Not keyed with code			

MISSION EQUIPMENT	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
RWR, APR-39(V)	Not installed			
IR CM, ALQ-144	Not installed			
Chaff dispenser, M-130	Not installed			
Cargo hook system	X			

HYDRAULIC CONTROL SYSTEM	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
Backup hydraulic pump	X			
Servo off 1st stage/PLT	X			
Servo off 2nd stage/PLT	X			
Servo off 1st stage/COPLT	X			
Servo off 2nd stage/COPLT	X			
Hydraulic leak test	X			
Tail servo	X			
Boost servos	X			

FUEL SYSTEM	No EMI Affect	EMI Affected Gnd Flt	Explanation
Fuel pump switch	X		
Fuel boost pump #1	X		
Fuel boost pump #2	X		
Fuel cont panel ESSS	X		
WARNING SYSTEM	No EMI Affect	EMI Affected Gnd Flt	Explanation
Low rotor RPM	X		
Master caution	X		
Caution advisory	X		
Fire warning	X		
AFCS	X		
Stabilator	X		
#1 engine out	X		
#2 engine out	X		
NAVIGATION INSTRUMENTS	No EMI Affect	EMI Affected Gnd Flt	Explanation
ADF	X		
Magnetic compass	X		
CONUS NAV, ARN-123	X		
Doppler, ASN-128	X		
Gyro mag compass (PLT)	X		
Gyro mag compass (COPLT)	X		
Compass cont panel, ASN-75	X		
HSI	X		
FLIGHT INSTRUMENTS	No EMI Affect	EMI Affected Gnd Flt	Explanation
Radar altimeter	X		
Stabilator pos indicator	X		
VSI	X		
CIS mode select	X		
SAS 1	X		
SAS 2	X		
FPS	X		
Trim	X		
Go-around enable	X		
Cyclic trim release	X		
Cyclic stick trim	X		
ALR encoder	X		

FLIGHT INSTRUMENTS (CONT)	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
HSI/VSI mode select (PLT)				
DPLR	X			
VOR/ILS	X			
BACK CRS	X			
FM HOME	X			
TURN RATE	X			
CRS HDG	X			
VERT GYRO	X			
BRG 2	X			
HSI/VSI Mode Select (COPLT)				
DPLR	X			
VOR/ILS	X			
BACK CRS	X			
FM HOME	X			
TURN RATE	X			
CRS HDG	X			
VERT GYRO	X			
BRG 2	X			
MISCELLANEOUS EQUIPMENT				
Blade deice	Not tested			Ambient tempera- ture was out of test lim- its.
Windshield anti-ice	X			
Pitot heat	X			
Vent blower	X			
Windshield wiper	X			
Heater	X			
APU	X			
Generator #1	X			
Generator #2	X			
Generator APU	X			
Air source heat start	X			
Tail wheel lock	X			
Gyro erect	X			

LIGHTING	No EMI Affect	EMI Affected		Explanation
		Gnd	Flt	
Cockpit utility	X			
Cockpit flood	X			
Cabin dome	X			
Search light	X			
Search light control	X			
Landing light	X			
Flt instr lights (PLT)	X			
Flt instr lights (COPLT)	X			
Nonflight instr lights	X			
Console lights, upper	X			
Console lights, lower	X			
Position lights	X			
Formation lights	X			
Anticollision lights	X			
NVG lighting	X			

**3.2.5 Battery life evaluation**

**Battery Life Evaluation  
Report Form**

**Nomenclature: Neonatal Monitor**  
**Manufacturer: Corometrics Medical Systems**  
**Model number: 506**  
**Serial number: AAC07901838**  
**Military item number: None**

**Options installed: None**

**Manufacturer battery life specification: Up to 4 hours.**

**Overall performance: Pass**

**Performance: Average 5 hours and 53 minutes operation on  
three trials from a fully charged battery.**

**Comments: Unit monitoring and ECG signal during test.**

### 3.2.6 Electrical safety test

#### Electrical Safety Test Report Form

Nomenclature: Neonatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: Corometrics Model 506  
Serial number: AAC07901838  
Military item number: None

Options installed: None

Date of test: 9 Jan 90

Performance:

Grounding conductor resistance (milliohms): 63

Leakage current - Case to ground (microamperes):

unit off, grounded, normal polarity	0.2
unit off, ungrounded, normal polarity	12.4
unit off, ungrounded, reverse polarity	18.8
unit on, grounded, normal polarity	0.1
unit on, ungrounded, normal polarity	12.5
unit on, ungrounded, reverse polarity	18.8

MAXIMUM LIMITS:

ground resistance (milliohms):	150
current (microamperes)	
current (grounded, type A unit):	10
current (ungrounded, type A unit):	100
current (grounded, type B unit):	50
current (ungrounded, type B unit):	500

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.7 Human factors evaluation

Human Factors Evaluation  
Report Form

Nomenclature: Neonatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: 506  
Serial number: AAC07901838  
Military item number: None

Options installed: None

Date of test: 9 Jan 90

Item configuration during test: Item prepared for operation.

Checklist for HFE

RESULTS

VISUAL DISPLAYS:

Satisfactory

display type, format, content  
location of displays  
indicator lights  
scalar displays  
color coding  
legends and labels  
cathode ray tubes  
counters  
flags, go-no-go, center-null indicators

Comments: Displays do not have brightness or intensity adjustments.

CONTROLS:

Unsatisfactory

location  
characteristics of controls  
labeling  
control - display relationships

Comments: Power switch is labeled "DISPLAY" not labeled "ON/OFF" or "POWER."

**TIME REQUIRED TO PREPARE FOR OPERATION (list in comment)**

**Comments: approximately 2 minutes.**

**MAINTAINABILITY:**

**Satisfactory**

component location  
component characteristics  
rests and stands  
covers, cases, access doors  
handles  
lubrication  
component mounting  
cord storage provisions  
external accessibility  
internal accessibility  
list special tools required  
list realistic inspection requirements  
list realistic inspection intervals

**Comments: No audible alarm for low battery condition.**

**CONDUCTORS:**

**Satisfactory**

binding and securing  
length  
protection  
routing  
conductor coding  
fabrication  
connectors

**Comments: None**

**FASTENERS:**

**Satisfactory**

access through inspection panel covers  
enclosure fasteners  
device mounting bolts and fasteners

**Comments: None**

**TEST POINTS:** **Satisfactory**

general  
location and mounting  
test point labeling and coding

Comments: None.

**TEST EQUIPMENT:** **Satisfactory**

general  
equipment self-test  
indicators (list in comments)  
controls  
positive indication of proper operation

Comments: Self test initiated when unit activated.

**FUSES AND CIRCUIT BREAKERS:** **Satisfactory**

external accessibility  
easy replacement or reset by operator

Comments: Fuses accessible on rear panel.

**LABELS AND CODING:** **Unsatisfactory**

placed above controls and displays  
near or on the items they identify  
not obscured by other equipment components  
describe the function of the items they identify  
readable from normal operating distance  
conspicuous placards adjacent to hazardous items

Comments: Power switch labeled "DISPLAY" and no illumination provided for controls.

**SAFETY:** **Satisfactory**

manual  
materials  
fire and explosive protection  
operator protection from mechanical hazards  
patient protection from mechanical hazards  
electrical safety (operator and patient)

Comments: No audible alarm for low battery condition.

### 3.2.8 Altitude test

#### Altitude Test Report Form

Nomenclature: Neonatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: Corometrics Model 506  
Serial number: AAC07901838  
Military item number: None

Options installed: None

Date of test: 11 Dec 89

Item configuration during test: Item sitting on chamber floor.

Performance test criteria: Proper display of simulated signals.

Ambient conditions outside chamber:

Temperature	70°F
Humidity	not recorded
Barometric pressure	1 atm

#### PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	ECG simulator
list connections to dummy loads	None
list unconnected terminals	Temp 2, BP, Computer I/O

#### IN-TEST DATA

Time of test start: 1335

#### POSTTEST DATA

Posttest performance check (complete check of item and accessories):

**Time of test end : 1450**

**Item functional (based on performance test criteria): Yes**

**Deviation from pretest : None**

**Comments on item setup or checks: None**

**Comments on test run (including interruptions): None**

**Comments on other data: None**

3.2.9 Vibration test

Vibration Test  
Report Form

Nomenclature: Neonatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: Corometrics Model 506  
Serial number: AAC07901838  
Military item number: None

Options installed: None

Date of test: 12 Dec 89

Item configuration during test: Item strapped down on vibration table fixture; ac operation.

Performance test criteria: Proper display of simulated signals.

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	ECG simulator
list connections to dummy loads	None
list unconnected terminals	Temp 2, BP, Computer I/O

Ambient conditions

Temperature	71°F
Humidity	61% RH
Barometric pressure	1 atm

IN-TEST DATA

Data and performance checks during test:

Time at first check:

X: 0940                      Y: 1245                      Z: 1220

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

**Time at second check:**

**X: 1020**

**Y: 1325**

**Z: 1300**

**Item functional (based on performance test criteria): Yes**

**Deviation from pretest: None**

**POSTTEST DATA**

**Time at test end:**

**X: 1030**

**Y: 1335**

**Z: 1310**

**Posttest performance check (complete check of item and accessories):**

**Item functional (based on performance test criteria): Yes**

**Item intact: Yes**

**Deviation from pretest: None**

**Comments on item setup or checks: None**

**Comments on test run (including interruptions): None**

**Comments on other data: Test times for the three axes are on different days.**

3.2.10 High temperature test

High Temperature Test  
(Equipment Operating)  
Report Form

Nomenclature: Neonatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: Corometrics Model 506  
Serial number: AAC07901838  
Military item number: None

Options installed: None

Date of test: 29 Aug 91

Item configuration during test: Unit was sitting on chamber floor.

Performance test criteria: Proper display of simulated signals.

Ambient conditions outside chamber:

Temperature	78°F
Humidity	50% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check :

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	ECG simulator
list connections to dummy loads	None
list unconnected terminals	Temp 2, BP, Computer I/O
distance from north wall (meters)	0.25
distance from south wall (meters)	0.41
distance from east wall (meters)	0.79
distance from west wall (meters)	1.75
distance from ceiling (meters)	1.20
distance from floor (meters)	0.7

IN-TEST DATA

Time of test start: 1130

Performance checks during test:

First check:

Time: 1200  
Temperature: 49°C  
Humidity: 15% RH  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria):  
Yes, all ok  
Deviation from pretest: None

**Second check:**

Time: 1230  
Temperature: 49°C  
Humidity: 15% RH  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria):  
Yes, all ok  
Deviation from pretest: None

**Third check:**

Time: 1300  
Temperature: 49°C  
Humidity: 15% RH  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria):  
Yes, all ok  
Deviation from pretest: None

**POSTTEST DATA**

**Posttest performance check:**  
(complete check of item and accessories)

Time of test end: 1330  
Item functional (based on performance test criteria):  
Yes, all ok  
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.11 High temperature storage test

High Temperature Test  
(Equipment in Storage)  
Report Form

Nomenclature: Neonatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: Corometrics Model 506  
Serial number: AAC07901838  
Military item number: None

Options installed: None

Date of test: 2 Jan 90

Item configuration during test: Sitting on chamber floor, in storage, not operating.

Performance test criteria: Proper display of simulated signals.

Ambient conditions outside chamber:

Temperature	24°C
Humidity	28% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	ECG simulator
list connections to dummy loads	None
list unconnected terminals	Temp 2, BP, Computer I/O
distance from north wall (meters)	0.25
distance from south wall (meters)	0.41
distance from east wall (meters)	0.79
distance from west wall (meters)	1.75
distance from ceiling (meters)	1.20
distance from floor (meters)	0.7

Time of test start: 0820

**POSTTEST DATA**

**Posttest performance check:**  
(complete check of item and accessories)

**Time of test end: 1153**  
**Item functional (based on performance test criteria): Yes**  
**Deviation from pretest: None**

**Comments on item setup or checks:**  
The unit was allowed to cool for 1 hour at ambient conditions before the posttest performance check was completed.

**Comments on test run (including interruptions): None**

**Comments on other data: None**

3.2.12 Low temperature test

Low Temperature Test  
(Equipment Operating)  
Report Form

Nomenclature: Neonatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: Corometrics Model 506  
Serial number: AAC07901838  
Military item number: None

Options installed: None

Date of test: 29 Dec 89

Item configuration during test: Sitting on chamber floor.

Performance test criteria: Proper display of simulated signals.

Ambient conditions outside chamber:

Temperature	23°C
Humidity	42% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Pass

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	ECG simulator
list connections to dummy loads	None
list unconnected terminals	Temp 2, BP, Computer I/O
distance from north wall (meters)	0.25
distance from south wall (meters)	0.41
distance from east wall (meters)	0.79
distance from west wall (meters)	1.75
distance from ceiling (meters)	1.20
distance from floor (meters)	0.7

Time of test start: 0815

Performance checks during test:

**First check:**

Time: 0845  
Temperature: 0°C  
Humidity: NA  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

**Second check:**

Time: 0915  
Temperature: 0°C  
Humidity: NA  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

**Third check:**

Time: 0945  
Temperature: 0°C  
Humidity: NA  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

**POSTTEST DATA**

Posttest performance check:  
(complete check of item and accessories)

Time of test end: 1030  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.13 Low temperature storage test

Low Temperature Test  
(Equipment in Storage)  
Report Form

Nomenclature: Neonatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: Corometrics Model 506  
Serial number: AAC07901838  
Military item number: None

Options installed: None

Date of test: 3 Jan 90

Item configuration during test: Sitting on chamber floor, not operating, in storage.

Performance test criteria: Proper display of simulated signals.

Ambient conditions outside chamber:

Temperature	24°C
Humidity	44% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	ECG simulator
list connections to dummy loads	None
list unconnected terminals	Temp 2, BP, Computer I/O
distance from north wall (meters)	0.25
distance from south wall (meters)	0.41
distance from east wall (meters)	0.79
distance from west wall (meters)	1.75
distance from ceiling (meters)	1.20
distance from floor (meters)	0.7

Time of test start: 0850  
Midtest time: 1230  
Midtest temperature: -46°C

**POSTTEST DATA**

**Posttest performance check:**  
(complete check of item and accessories)

Time of test end: 1600  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.14 Humidity test

Humidity Test  
Report Form

Nomenclature: Neohatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: Model 508  
Serial number: AAC07901838  
Military item number: None

Options installed: None

Date of test: 29 Dec 89

Item configuration during test: The unit was sitting on the chamber floor, operating on ac power.

Performance test criteria: Proper display of simulated signals.

Ambient conditions outside chamber:

Temperature	25°C
Humidity	58% RH
Barometric pressure	1 atm

PRETEST DATA

Pretest performance check:

Item functional (based on performance test criteria): Yes

Installation of item in test facility:

list connections to power	120 Vac
list connections to simulators	ECG simulator
list connections to dummy loads	None
list unconnected terminals	Temp 2, BP, Computer I/O
distance from north wall (meters)	0.25
distance from south wall (meters)	0.41
distance from east wall (meters)	0.79
distance from west wall (meters)	1.75
distance from ceiling (meters)	1.20
distance from floor (meters)	0.7

IN-TEST DATA

Time of test start: 1120

**Performance checks during test:**

**First check:**

Time: 1230  
Temperature: 29.5°C  
Humidity: 95% RH  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

**Second check:**

Time: 1315  
Temperature: 29.5°C  
Humidity: 95% RH  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

**Third check:**

Time: 1400  
Temperature: 29.5°C  
Humidity: 95% RH  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

**Fourth check:**

Time: 1445  
Temperature: 29.5°C  
Humidity: 95% RH  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

**Fifth check:**

Time: 1520  
Temperature: 29.5°C  
Humidity: 95% RH  
Barometric pressure: 1 atm  
Item functional (based on performance test criteria): Yes  
Deviation from pretest: None

**POSTTEST DATA**

**Posttest performance check:**

(complete check of item and accessories)

Time of test end: 1610

Item functional (based on performance test criteria): Yes

Deviation from pretest: None

Comments on item setup or checks: None

Comments on test run (including interruptions): None

Comments on other data: None

3.2.15 Electromagnetic characteristics test

\*\*\*\*\*

Electromagnetic Characteristics Testing  
Evaluation of Performance

\*\*\*\*\*

T & E Item Number: 19

Date: 17 Jan 90

Nomenclature: Neonatal Monitor  
Manufacturer: Corometrics Medical Systems  
Model number: 506  
Serial number: AAC07901838  
Military item number: NA

\*\*\*\*\*

Conducted Emissions Tests

- CE01      Testing configuration(s): NA  
            Performance (pass/fail): NA  
  
            Comments: No dc conductors
- CE02      Testing configuration(s): Operating on copper  
            work bench.  
            Performance (pass/fail): Pass  
  
            Comments: No signal failures.
- CE04      Testing configuration(s): Operating on copper  
            work bench.  
            Performance (pass/fail): Pass  
  
            Comments: No signal failures.

Conducted Susceptibility Tests

- CS02      Testing configuration(s): Operating on test  
            bench, connected to test jig.  
            Performance (pass/fail): Pass  
  
            Comments: Not susceptible to test signals on  
                    power conductors.

CS06      Testing configuration(s): Operating on counter top.  
 Performance (pass/fail): Pass  
 Comments: Not susceptible to test spikes

**Radiated Emissions Tests**

RE02      Testing configuration(s): Operating on wooden test stand in the EMC chamber, ac power.  
 Performance (pass/fail): Fail  
 Comments: Emission data:

	<u>Frequency range</u>	<u>Emission exceeding standard</u>
ac:	0.066 - 107.2 MHz	0.2 - 25.5 dB (NB)
	0.014 - 0.15 MHz	3.2 - 10.4 dB (BB)
battery:		
	0.015 - 107.275 MHz	0.4 - 36.5 dB (NB)
	0.014 - 0.55 MHz	4.3 - 25.1 dB (BB)
	1.025 MHz	1.3 dB (BB)
	1.672 MHz	8.5 dB (BB)

**Radiated Susceptibility Tests**

RS03      Testing configuration(s): Operating on the wooden test stand in the EMC chamber.  
 Performance (pass/fail): Fail  
 Comments: Susceptibility data:

	<u>Frequency</u>	<u>Field strength</u>
ac:	6.7 - 7.7 MHz	0.53 - 1.67 V/m
	14.8 - 22.4 MHz	0.28 - 3.34 V/m
	30.0 - 176.2 MHz	0.63 - 7.94 V/m
battery:		
	20.4 - 24.8 MHz	0.35 - 4.21 V/m
	30.0 - 183 MHz	0.53 - 8.91 V/m

### 3.3 CRITERIA, SIGNIFICANT PROBLEMS, AND SUGGESTED IMPROVEMENTS

#### 3.3.1 Criteria

<u>Item</u>			<u>Applicable</u>
<u>No.</u>	<u>Criteria (source)</u>	<u>Remarks</u>	<u>subparagraph</u>
1	The physical inventory is conducted solely for investigation and documentation.	NA	2.1.2.1
2	The Corometrics Model 506 will display consistent and accurate performance.	met	2.1.2.2
3	Verify manufacturer's specified full power internal battery life expectancy of 4 hours.	met	2.2.2
4	The Corometrics Model 506 will meet the limits established in NFPA 99 for electrical safety of medical equipment.	met	2.3.2
5	The Corometrics Model 506 will be rated satisfactory in all major categories of the evaluation. These include: Visual displays, controls, maintainability, conductors, fasteners, test points, test equipment, fuses and circuit breakers, labels and coding, and safety.	partially met	2.4.2
6	The Corometrics Model 506 will demonstrate proper operation - while exposed to an altitude equivalency of 15,000 feet above sea level.	met	2.5.2
7	The Corometrics Model 506 will remain operational while exposed to vibrational stresses.	met	2.6.2
8	The Corometrics Model 506 will remain operational during the high temperature operation check.	met	2.7.2.1

9	The Corometrics Model 506 will remain operational after the high temperature storage.	met	2.7.2.2
10	The Corometrics Model 506 will remain operational during the low temperature operation check.	met	2.8.2.1
11	The Corometrics Model 506 will remain operational after the low temperature storage.	met	2.8.2.2
12	The Corometrics Model 506 will remain operational while exposed to a high humidity.	met	2.9.2
13	The Corometrics Model 506 will not produce emissions in excess of the limits set forth in MIL-STD-461A Notice 4, paragraph 6.13.	partially met	2.10.2.1
14	The Corometrics Model 506 will not malfunction when it is subjected to radiated fields as specified in MIL-STD-461A, Notice 4, paragraph 6.20.	partially met	2.10.2.2
15	The Corometrics Model 506 will not conduct emissions in excess of the limits set forth in MIL-STD-461A, Notice 4, paragraph 6.2.	met	2.10.2.3
16	The Corometrics Model 506 will not malfunction when it is subjected to conducted emissions as specified in MIL-STD-461A, Notice 4, paragraphs 6.7 and 6.10.	met	2.10.2.4
17	The flight surgeon will be able to operate the Corometrics Model 506 without physical or functional restrictions aboard the aircraft.	met	2.11.2.1

- |    |   |     |          |
|----|---|-----|----------|
| 18 | The Corometrics Model 506 will not radiate EMI to disrupt or interfere with the other equipment or systems aboard the aircraft. | met | 2.12.2.2 |
| 19 | The aircraft will not radiate EMI to disrupt or interfere with the Corometrics Model 506.                                       | met | 2.12.2.3 |

3.3.2 Significant problems which require corrective action

None

3.3.3 Suggested improvements

None

### 3.4 REFERENCES

3.4.1 Department of Defense. 1971. EMI characteristics, requirements for equipment. Washington, DC. MIL-STD-461A, Notice 4. February.

3.4.2 Department of Defense. 1971. EMI characteristics, measurement of. Washington, DC. MIL-STD-462, Notice 3. February.

3.4.3 Department of Defense. 1983. Environmental test methods and engineering guidelines. Washington, DC. MIL-STD-810D. July.

3.4.4 Department of the Army. 1987. Maintenance management procedures for medical equipment. Washington, DC. TB 38-750-2. April.

3.4.5 Underwriters Laboratory's, Inc. 1978. Standard for safety, medical and dental equipment. Chicago, Illinois. UL-544.

3.4.6 Department of Defense. 1989. Human engineering design criteria for military systems, equipment, and facilities. Washington, DC. MIL-STD-1472D. March.

3.4.7 Association for the Advancement of Medical Instruments. 1988. Human factors engineering guidelines and preferred practices for the design of medical devices. Arlington, Virginia. AAMI-HE-1988. February.

3.4.8 National Fire Protection Association. 1987. Standard for health care facilities. Quincy, Massachusetts. NFPA 99. February.

3.4.9 Department of the Army. 1982. Environmental protection and enhancement. Washington, DC. AR 200-1. June.

### 3.5 ABBREVIATIONS

ac	alternate current
AVSCOM	Army Aviation Systems Command
AWR	airworthiness release
BB	broadband
CAAF	Cairns Army Airfield
dc	direct current
EMC	electromagnetic compatibility
EMI	electromagnetic interference
fpm	feet per minute
GFE	government furnished equipment
Gpk	gravity, peak
G(rms)	gravity (root mean square)
Hz	hertz
IAW	in accordance with
ITOP	in-flight test operating procedure
IV	intravenous
kHz	kilohertz
LCD	liquid crystal display
LED	light emitting diode
LISN	line impedance stabilization network
MEDEVAC	medical evacuation
MHz	megahertz
MIL-STD	military standard
mL	milliliter
mm	millimeter
mmHg	millimeters of Mercury
MSL	mean sea level
NFPA	National Fire Prevention Association
NB	narrowband
NBC	nuclear, biological and chemical
NOE	nap-of-the-earth
NVG	night vision goggle
RF	radio frequency
RFI	radio frequency interference
RH	relative humidity

TB  
TFT  
T & E

technical bulletin  
technical feasibility testing  
test and evaluation

UES  
USAARL

Universal Energy Systems, Inc.  
U.S. Army Aeromedical Research Laboratory

V/m

volts per meter

### 3.6 LIST OF MANUFACTURERS

- 3.6.1 Corometrics Medical Systems, Inc.  
61 Barnes Park Road North  
Wallingford, CN 06492
- 3.6.2 Neurodyne-Dempsey, Inc.  
200 Arrowhead Drive  
Carson City, NV 89701
- 3.6.3 Tenney Engineering, Inc.  
1090 Springfield Road  
P.O. box 3142  
Union, NJ 07083
- 3.6.4 Unholtz-Dickey Corporation  
6 Brookside Drive  
Wallingford, CT 06492
- 3.6.5 Solar Electronics Company  
901 North Highland Avenue  
Hollywood, CA 90038
- 3.6.6 Tektronix, Inc.  
P.O. Box 500  
Beaverton, OR 97077